# MOVIVE™

## INSTRUCTIONS FOR USE

...▶ Take the next step by reviewing this guide  
and learn how Movive can help


Rx Only

# Table of Contents

---

<b>What Is Move?</b> .....	<b>2</b>
How Does it Work?.....	4
Move Components.....	5
Where Should I Walk? .....	8
<b>First Time Set-Up</b> .....	<b>9</b>
How to Charge Move .....	11
Turn on Touchscreen Device and Follow Instructions.....	13
Attach the Foot Sensors and Turn on the Headset .....	16
<b>Safety Information</b> .....	<b>19</b>
Warnings .....	20
Precautions .....	22
<b>Move Session From Start to Finish</b> .....	<b>24</b>
How to Start a Session.....	25
How to Confirm the System Is Connected.....	26
How to Select Your Music .....	28
How your Foot Sensors Calibrate .....	30
Walking a Session .....	31
How to End a Session .....	33
How to Store & Charge Move .....	35
<b>Troubleshooting and Support</b> .....	<b>36</b>
<b>Spanish Instructions for Use</b> .....	<b>37</b>

# What Is Moveive?

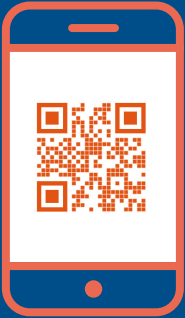Watch a welcome video by visiting [Moveive.com/welcome](https://Moveive.com/welcome) or by scanning this QR code with your personal smartphone.

Additionally, you can watch it anytime directly on the Moveive touchscreen device.

## **INDICATIONS FOR USE:**

Moveive is intended to support gait rehabilitation and motor-function for people living with Parkinson's disease.

# WELCOME!

You are about to take the **next step in your rehabilitation journey** and Move is at your side.


Move uses a combination of music, foot sensors, and advanced software to support walking rehabilitation. The intervention is based on the latest neuroscience research shown to support physical rehabilitation of adults with walking impairments from Parkinson's Disease.

## You Will:

---

Clip the 2 foot sensors to your shoes.

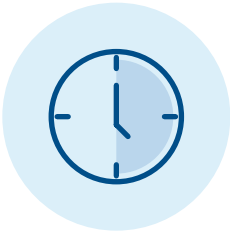Use the touchscreen device to guide you through a 30-minute walking session.

Walk a session.

...▶ Please use Move as prescribed by your doctor.

## How Does Move Work?

---

Move automatically adjusts music to improve walking.

- The system uses **2 foot sensors** to measure the way you walk, and a **touchscreen** device to guide each session.
- The system uses the information from the foot sensors to set the starting speed of the music.
- **You will hear audio** while you walk.


**During each session, the music speed will slowly increase or decrease based on how you are walking.**

- You can use Move in your home.
- To stay up to date, Move periodically needs access to WiFi or a cellular data connection. The touchscreen device comes with a preloaded, cellular plan. (See page 14 for details)


# Movive Rehabilitation System Components

---


**1x Headset**

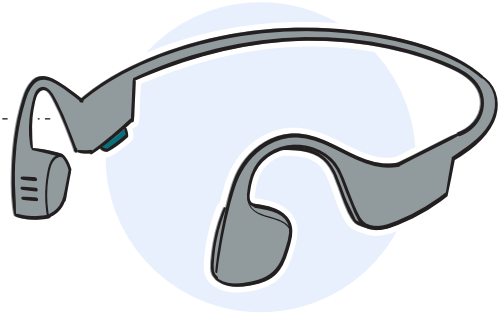**2x Foot sensors  
(left and right)**

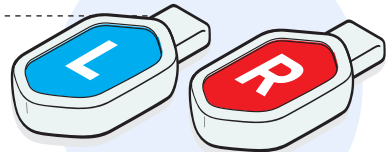**1x Touchscreen  
device**

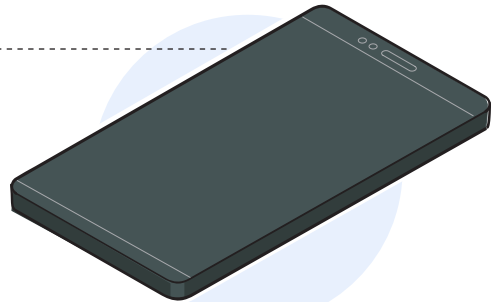
**1x Charging  
equipment**

1x Hub

3x Cables

1x AC Adapter and Cable

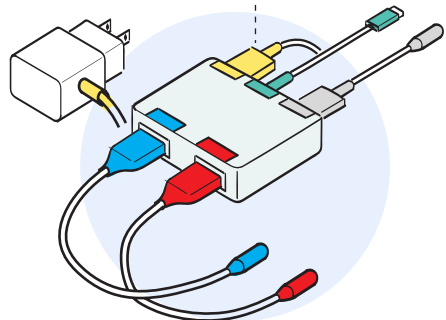\*The hardware in your kit may slightly differ from what is shown


## 1 Left and Right Foot Sensors

- By wearing the foot sensors, Movive is able to evaluate how you walk.
- The foot sensors 'turn on' just by picking them up. You don't have to press a button.
- The foot sensors are wireless. They are already connected (via Bluetooth) to the touchscreen device. They do not require any additional setup. **Please charge them after every use.**

## 2 Touchscreen Device


- The touchscreen device can only be used for Movive. The software is already installed on the device. You do not need to download anything.
- The touchscreen device will show you how to start, pause, and end a session.
- Keep the touchscreen device within arm's length during each session (10 feet) to maintain the connection with the sensors.

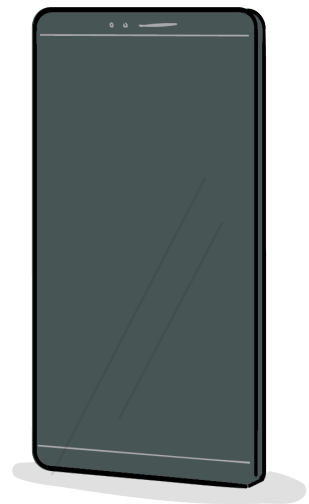\*You can listen to the audio from the speakers in the touchscreen device.


### 3 Headset

- The headset can be used to listen to the music while you walk.
- The wireless headset is already connected (via Bluetooth) to the touchscreen device. **You only have to turn it on.**
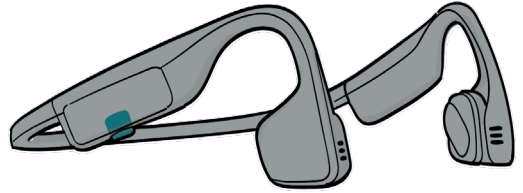
### 4 Charging Equipment

- To charge the foot sensors, touchscreen device, and headset, you must use the Move charging equipment.
- The data from your session will upload to secure cloud storage while you charge the system.
- The system is designed to charge and upload data **only with this charging equipment.** Switching out cables may impact system performance.
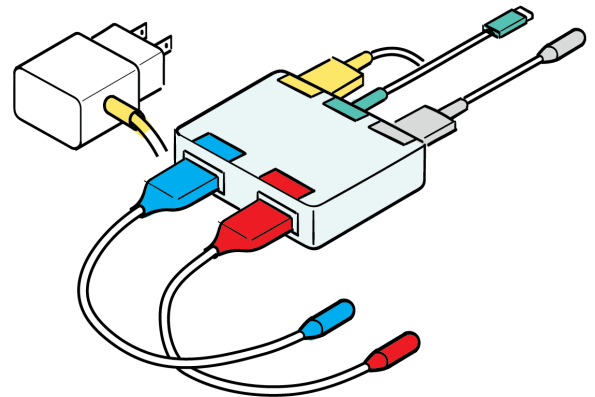
## Where Should I Walk?

---

Use Move in locations where the foot sensors can most accurately detect steps:

Flat and level surfaces


Locations where you can walk at least 20 steps in a straight line before turning. This enables the foot sensors to get a good read of your steps.

**DO NOT** use Move on a treadmill.

...▶ Please see page 19 section for full details on safety information.

# First-time Set Up

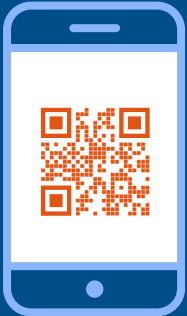Watch a welcome video by visiting [Movive.com/welcome](https://Movive.com/welcome) or by scanning this QR code with your personal smartphone.

Additionally, you can watch it anytime directly on the Movive touchscreen device.

## Set up Move in 3 Simple Steps:

---

Choose a safe place to walk that's flat and level.

Turn on the touchscreen device and attach the foot sensors.


When the music starts, all you need to do is walk!

...▶ See page 25 for details on how to start a session.

# Charge the System

Before you begin, **charge the system for at least 1 hour.** Ensure all charging cables are securely connected together as pictured below.

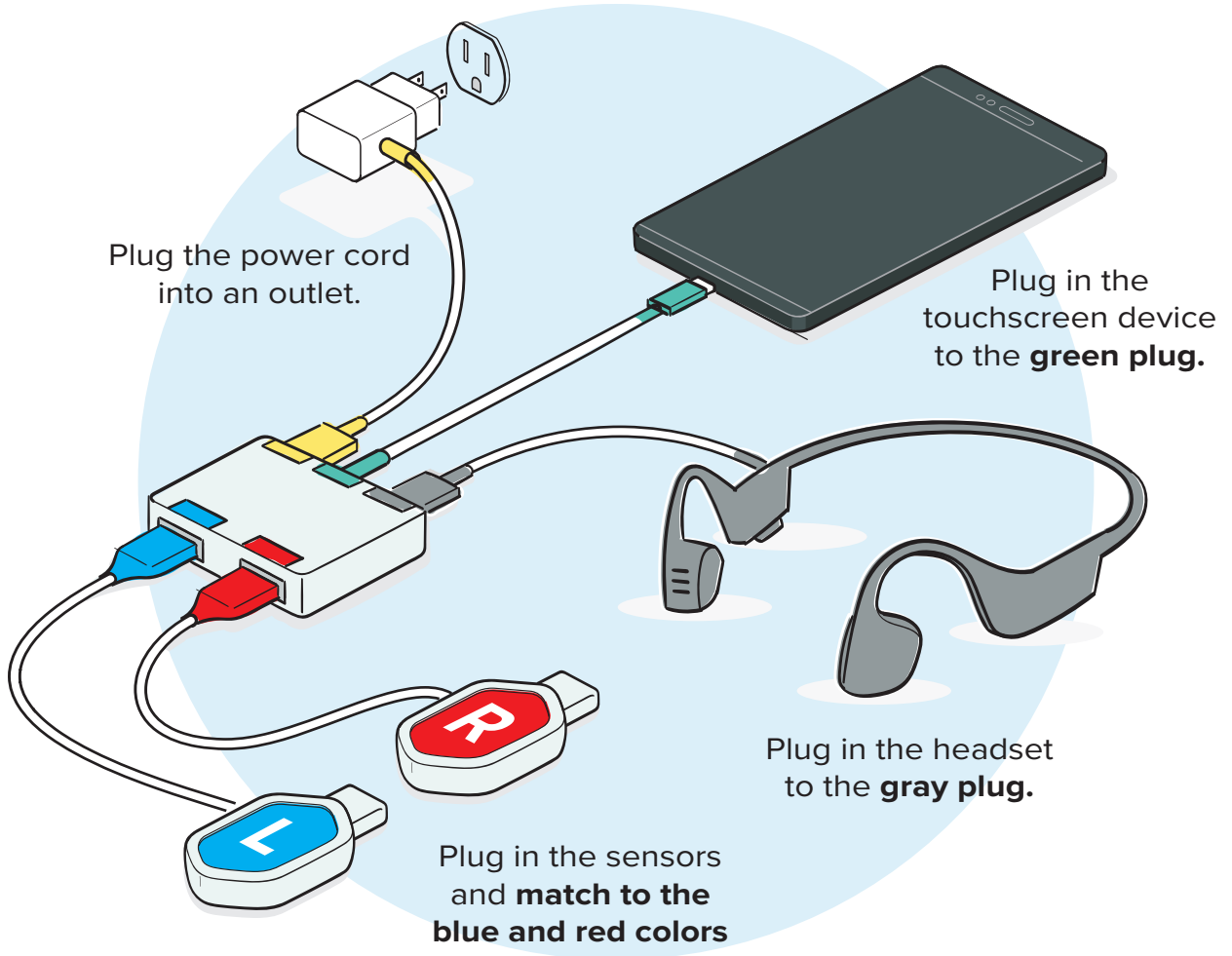
Always charge before use.


The cables are magnetic to help each piece stay in place to charge.

## Check Battery Levels

- When plugged into the charging equipment, the touchscreen device will display the battery status.

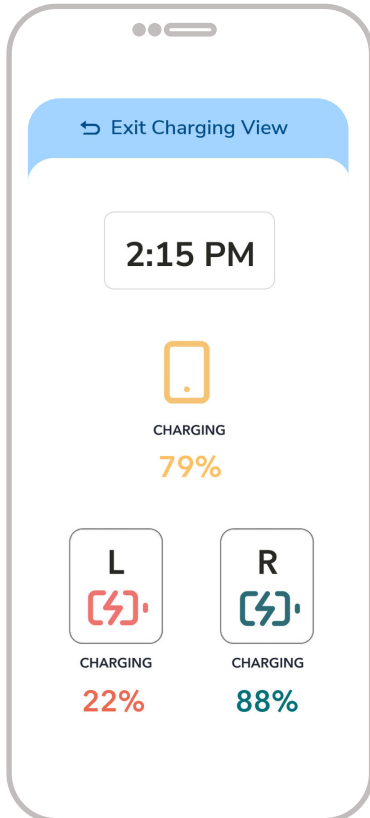**Charging in progress**

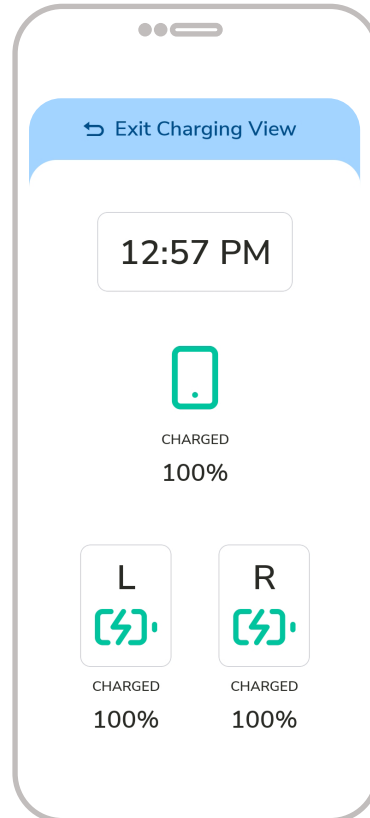
**System fully charged**

- Before each session, the touchscreen device **automatically checks all battery levels** before you begin walking.
  - » If the battery levels are sufficient, you will see a confirmation that the system is ready to use.
  - » If the battery levels are too low to start a session, you will be asked to charge the system.

# Turn on the Touchscreen Device and Follow the Instructions

---

## Power on the Touchscreen Device

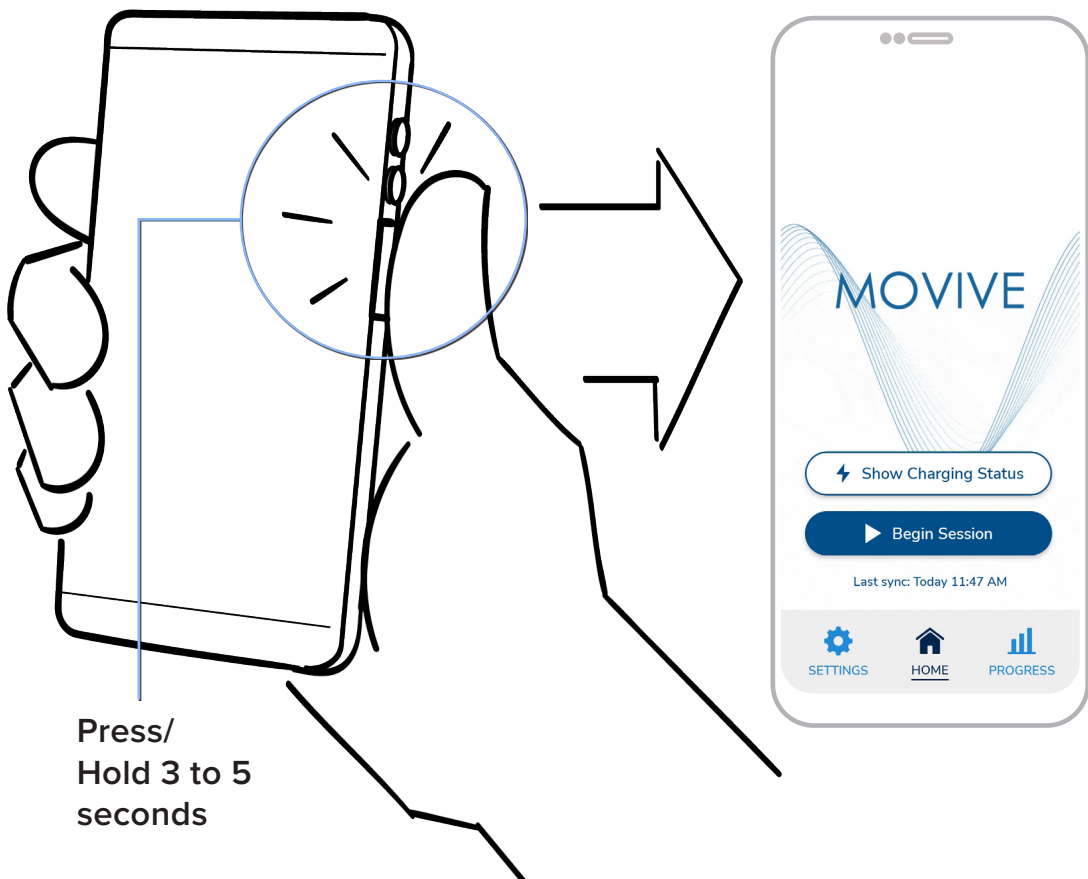
- Turn on by pressing and holding the power button until you feel a vibration, about 3 to 5 seconds.
- Press the  button.

**Moveive relies on a good WiFi or cellular data connection to get updates, new music playlists, and to save your walking data.**


## Choose a Data Connection

- During first-time setup, you will be asked if you would like to **set up WiFi**. If you have access to WiFi, it's recommended to use that for your connection.
- The touchscreen device already has a **cellular plan ready to use**. This means it usually has a good data connection without extra set up.
- If you'd like to connect to a WiFi network later, please go to the **“Settings”** menu on the home screen.”
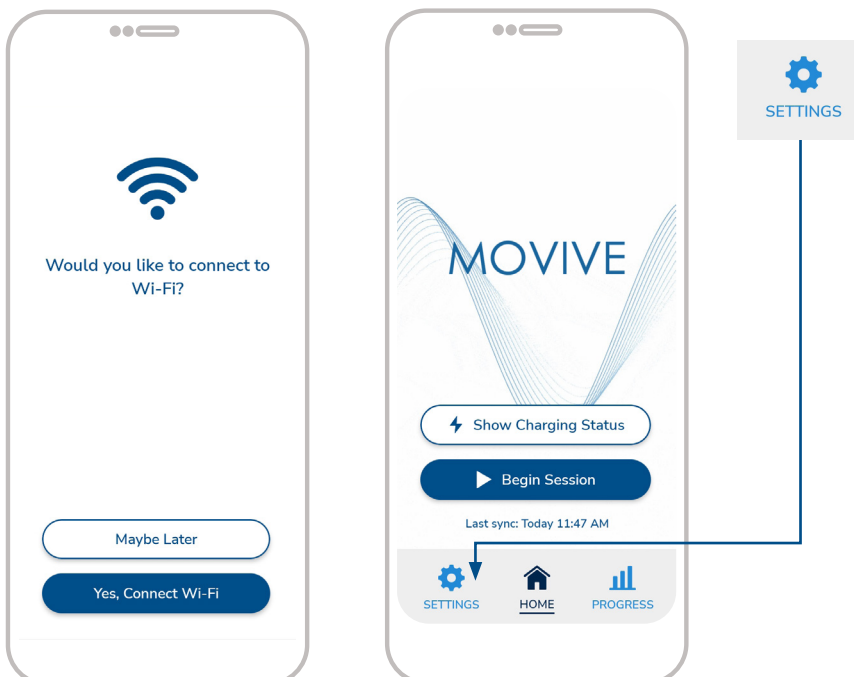
## View Setup Reminders


Before you begin a session, the touchscreen device will guide you through setup reminders. You will be given the option to **Watch a video** about how the system works.
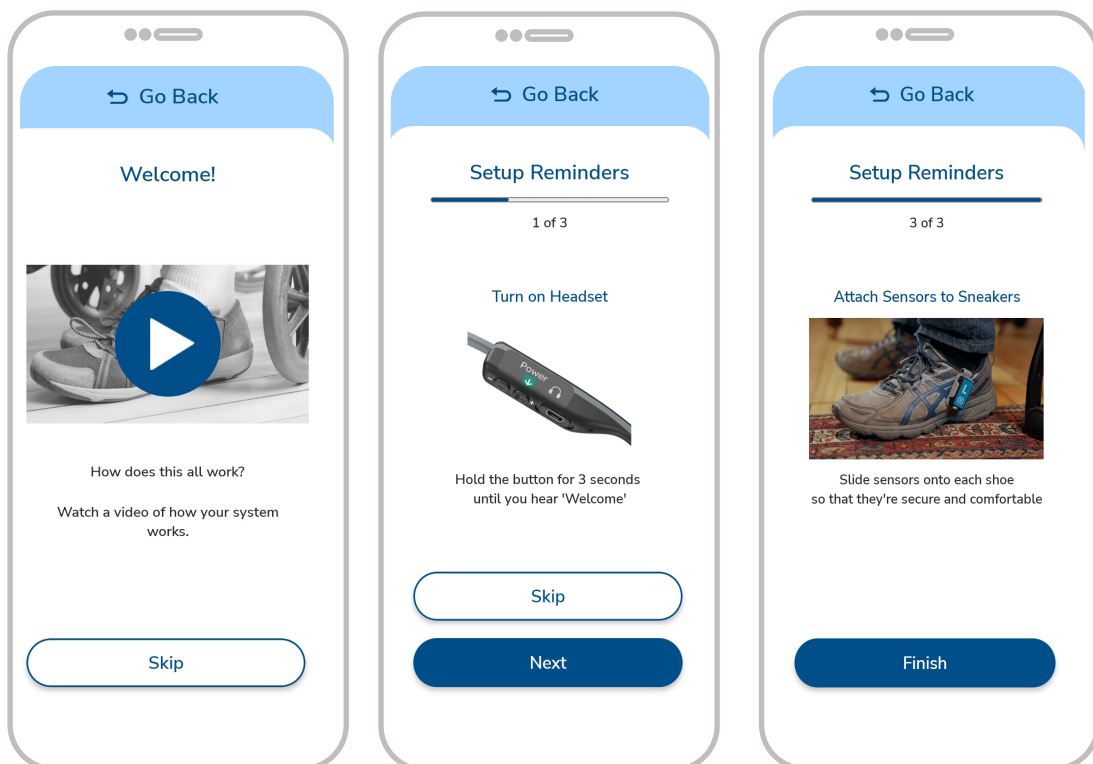
# Attach the Foot Sensors and Turn on the Headset

## Slide Foot Sensors Onto Your Shoes

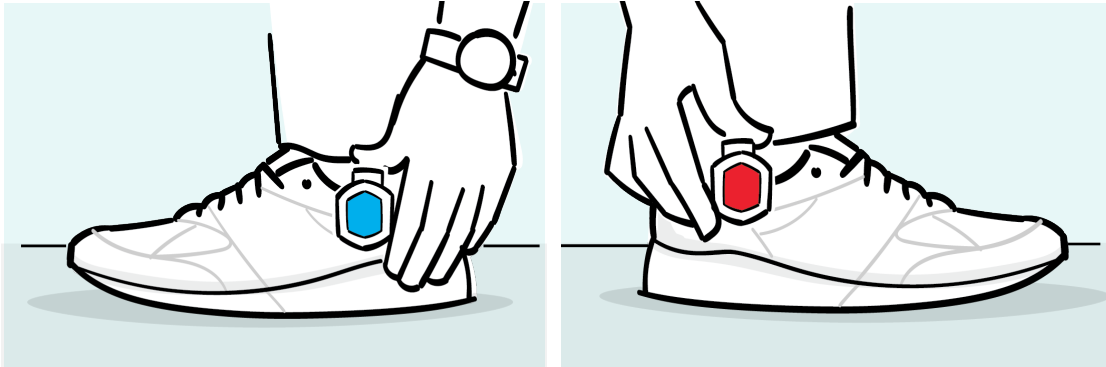
Slide each sensor onto the appropriate shoe.

- Foot sensors are labeled “L” and “R” for **left** and **right** feet, accordingly. Having a dedicated sensor for each foot helps the system better evaluate your walking.
- Foot sensors can go **anywhere on your shoes** as long as they feel comfortable and secure. It’s recommended to put them on the outside face of your shoes.
- Try standing with the foot sensors on to **make sure they feel comfortable** and won’t fall off while you walk.
- When charged, the **foot sensors turn on just by picking them up. You don’t have to press a button to turn them on.**

**Put on foot sensors from a stable position, such as sitting down, or ask someone for help to prevent loss of balance and risk of a fall.**

## Verify the Foot Sensors Have Connected

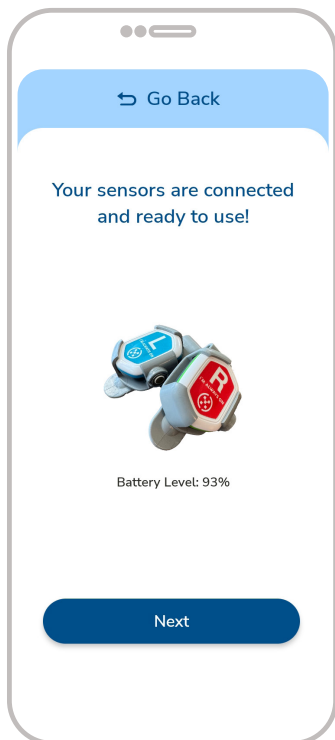
- The foot sensors are connected to the touchscreen device wirelessly (via Bluetooth). **They are already connected** and do not need additional set up.
- Once you start a session, the **touchscreen device will confirm the foot sensors are connected and charged** before you walk.

Keep the touchscreen device close to the foot sensors (within 10 feet) during your session to maintain the connection.

## Turn on the Headset

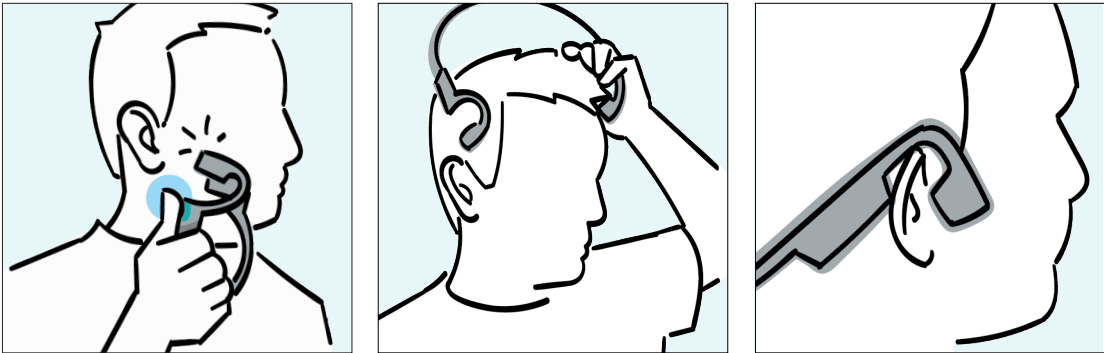- Turn on the headset by holding down the **small button near the green arrow**  for 3 to 5 seconds.
- While you press down the button, **hold the headset close to your ear**. You will hear a **confirmation sound** that it turned on.
- When the headset is turned on, a **blue indicator light will briefly light up**.
- The headset **speakers should rest comfortably in front of your ear**. They are designed to rest on your cheekbone rather than inside of your ear. This will help you to be able to hear your environment while you walk.


# Safety Information

This section lists general warnings and precautions related to using Movive.

## **WARNING**

Describes **potential serious** adverse reactions or safety hazards.
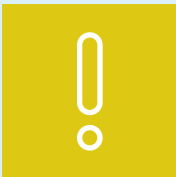
## **PRECAUTION**

Describes potentially hazardous situations that may result in minor to moderate injury or damage to the equipment.

## Warnings

---

### To prevent risk of injury from a slip, trip, fall, or collision:

- Attach and remove the foot sensors from your shoes from a stable position, such as sitting down.
- Use Move in a safe walking area.
  - » **DO NOT** use in poorly-lit areas.
  - » **DO NOT** use in areas with obstacles or potential hazards in your walking path, such as pets, pedestrians, or traffic.
  - » **DO NOT** use in areas with uneven, wet, or icy terrain.
  - » **DO NOT** use with the volume very loud or in loud surroundings because it may hide safety cues, such as emergency sirens or cars honking.
  - » **DO NOT** walk while looking at the touchscreen device.
  - » **DO NOT** continue using Move if walking becomes too difficult. Take breaks as needed. A session can be paused and resumed using the touchscreen device at any time throughout the session.
- If you ever become uncomfortable or feel unsafe with the speed of the music, slow down. The speed of the music will adjust back to the starting pace of your walking.

## Warnings *continued*

---

- If you experience muscle or joint pain or discomfort that you normally do not have, stop using Moveive and consult your doctor.
- If you experience any abnormal symptoms such as nausea, shortness of breath or excessive fatigue, stop using Moveive and consult your doctor.
- **DO NOT** use Moveive if any of the system components are damaged (to prevent risk of electric shock).
- **DO NOT** try to repair your Moveive system or charging equipment. Contact Moveive Support if you experience a technical problem not covered in this guide.

## Precautions

---

### Use Move only as directed by these Instructions for Use and your doctor.

- Wear footwear that is well fitted and supportive.
- Wear calf-height socks to prevent possible discomfort from the sensor clips against your skin.
- Check that the foot sensors are securely attached to your shoes. If a sensor falls off during a session, retrieve it if you can do so safely.

### To help Move work properly and ensure accurate step detection from the foot sensors:

- Use Move in areas where you can walk at least 20 steps in a straight line before turning.
- **AVOID** locations with high noise levels that may limit your ability to hear system audio.
- **AVOID** uneven walking surfaces such as frequent stairs, ramps, or hiking trails.
- **DO NOT** use Move on a treadmill.
- Move should only be used by the patient for whom it was prescribed. It is recommended that the device be safely stored to reduce the risk of unauthorized access.

## Precautions *continued*


---

### To prevent damage and/or performance issues with Moveive:

- Charge the foot sensors, touchscreen device, and headset for at least an hour before the start of each session.
- Only charge Moveive with the provided charging equipment. The components provided with Moveive use USB-style 5V/1A power blocks that must be plugged into 120V/60Hz power outlets.
- **DO NOT** use Moveive when it is charging.
- Use a bag to carry the Moveive system and charging equipment.
- Store in Ambient Conditions. **DO NOT** store in extremely hot, cold, humid, wet, or bright conditions to ensure proper system performance.

# Movive Session From Start to Finish

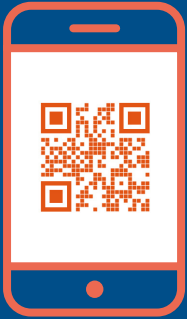Watch a welcome video by visiting [Movive.com/welcome](https://www.movive.com/welcome) or by scanning this QR code with your personal smartphone.

Additionally, you can watch it anytime directly on the Movive touchscreen device.

To start a walking session with Movive, put on the foot sensors and headset, and power on the touchscreen device. Next you will:

**Press the 'Begin Session' Button**

---

**Confirm the System Is Connected**

---

**Walk a Session**


---

**Charge the System**

# How to Start a Session

---

Start a session at any time from the touchscreen device's home screen.

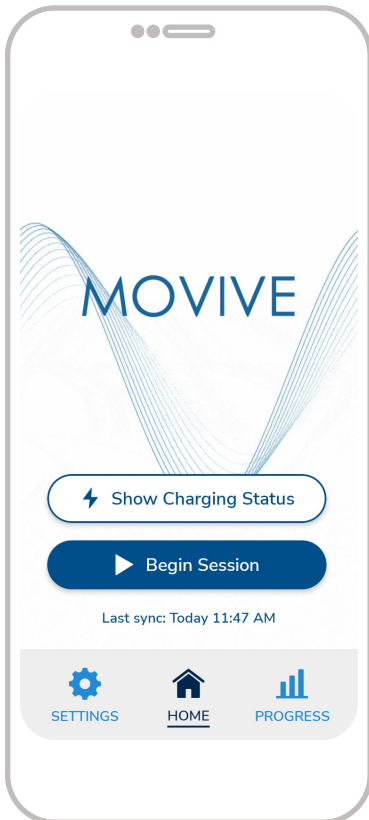Press the


button to

begin.


# How to Confirm the System Is Connected

---


## Check the Connection

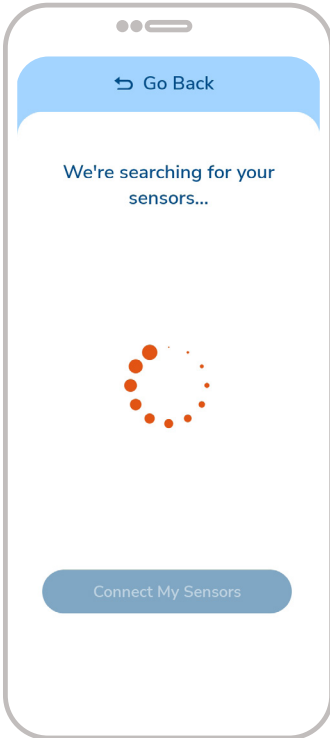The touchscreen device will automatically check whether the foot sensors and headset are connected and charged.

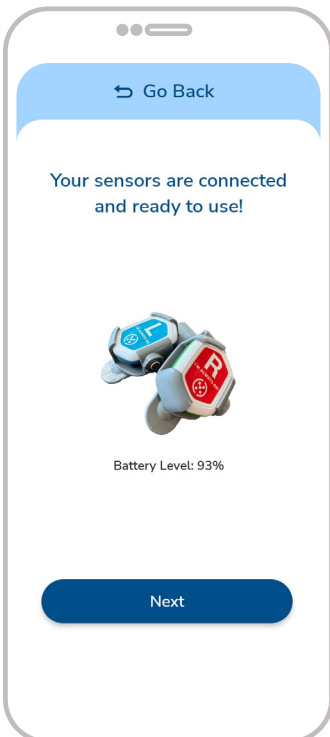
If the foot sensors and headset are ready to use, you will see a confirmation.

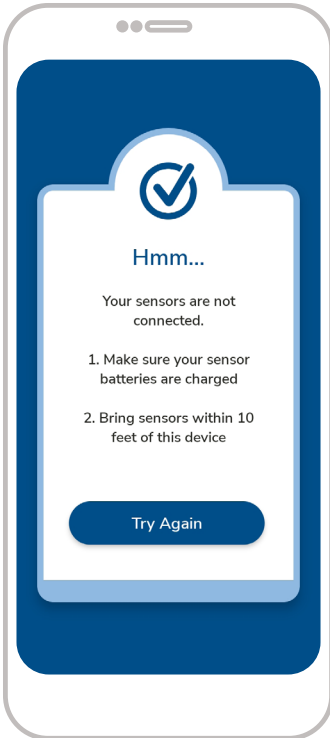An issue connecting could be caused by the foot sensors not having enough charge. An issue connecting could also be caused by the foot sensors being more than 10 feet from the touchscreen device.

**If the issue connecting foot sensors continues, please try restarting the touchscreen device by holding the power button for 3-5 seconds.**

## Check the Volume

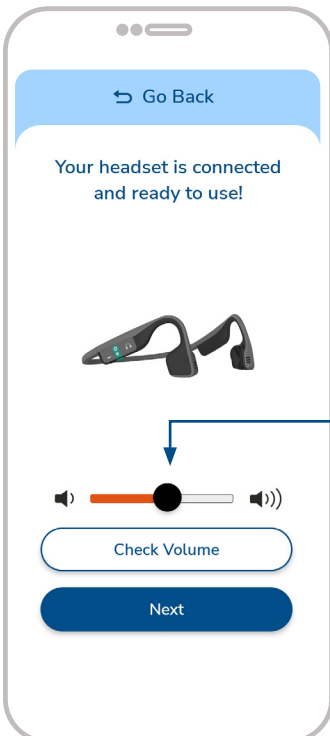
Test the volume by pressing the **Check Volume** button.

This will play a short music sample.

**Adjust the volume** to a comfortable level by using either the sliding scale on the touchscreen or the volume buttons on the side of the device.


# How to Select Your Music

---

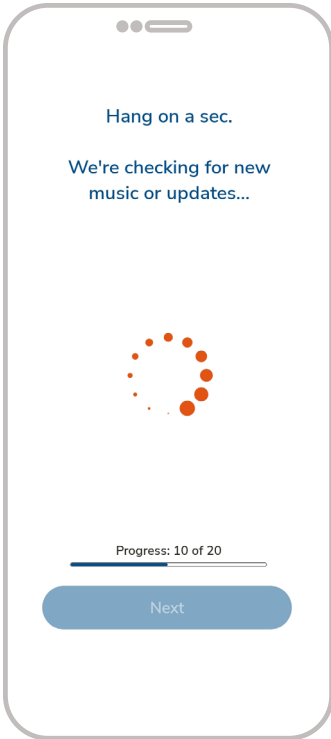
There is music preloaded on the touchscreen device, **ready to use out of the box.**

Movive will automatically check for new music or updates.

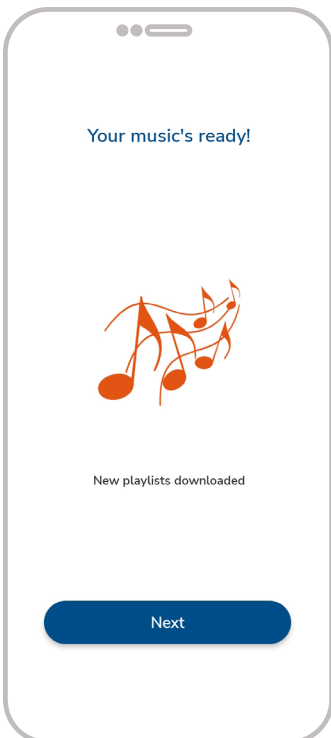If new music was available, you will receive the following message: **New playlists downloaded.**

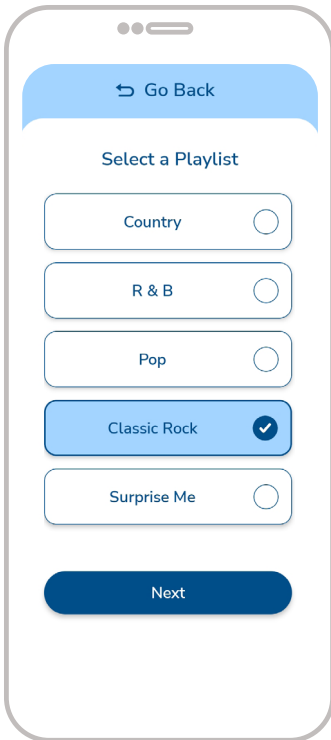Before each session, you will be prompted to **select your preferred playlist**.

You can select a different playlist for each session, but playlists **cannot be changed during a session**.

## Confirm You're Ready to Walk

This confirmation screen shows when the system is ready.

- Press the **I'm Ready to Walk** button
- It's recommended that you **put the touchscreen device away in a pocket** during each session.
- Begin walking at a **comfortable pace** so that the foot sensors can establish your starting pace.

# How Your Foot Sensors Calibrate

---

## Sensor Calibration

- **Each session begins by asking you to walk at a comfortable pace for about 60 seconds.** During this time, Move learns about your walking and personalizes the session. This is called “sensor calibration.”
- **Sensor calibration helps Move determine what speed to start your music.** Every day your walking pace can be different, so Move always starts with a new calibration.
- **It’s important to start your walk at a pace that you’d be comfortable walking at for 30 minutes.** That way, if the session ever becomes too difficult, Move can return to this starting pace. During the session, you can always take breaks as needed.

## Walking a Session

---

When the music starts to play, all you need to do is walk.

- Over the course of a 30-minute session, the music will gradually change speed and become more challenging. **Continue to listen to the music** and simply walk.
- Remember, the goal of the system is to **gradually change your speed**, and at times that may gently push you past your comfort zone.
- Moveive is also monitoring whether a pace seems too difficult based on the quality of your walking and **will respond as needed. If you ever become uncomfortable or feel unsafe with the speed of the music, slow down.** The speed of the music will adjust back to the starting pace of your walking.


The goal of the therapy is to challenge you past your comfort zone!


## Take Breaks, if Needed

- **It is okay to take a break for a few minutes** if you feel tired or need to make some adjustments to continue walking more comfortably.
- The goal is to **walk as much as possible** within the 30-minute session.

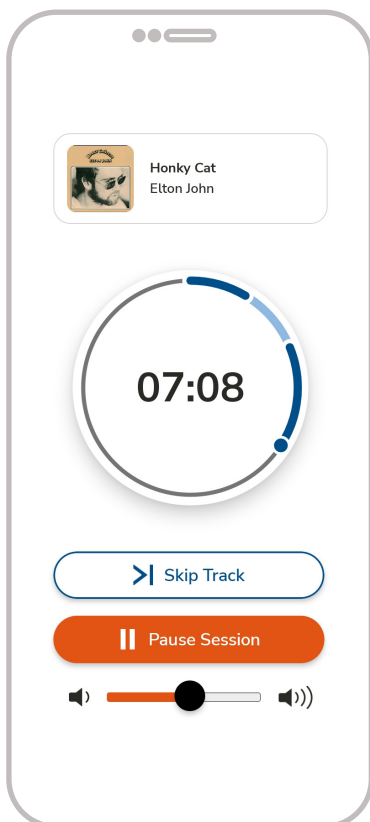If you would like the music to stop playing during a break, simply press 

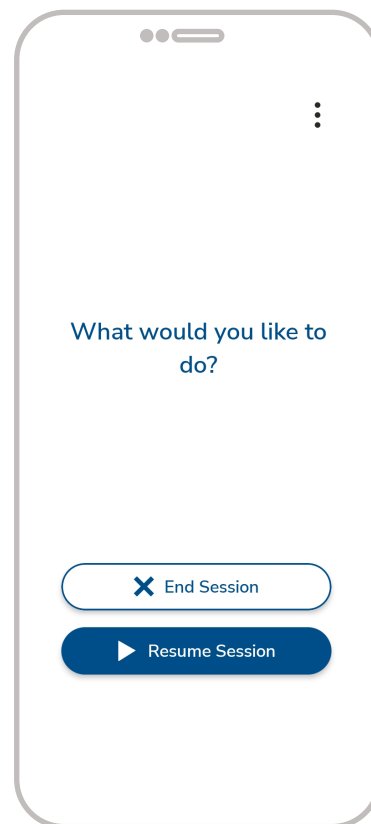
Press  to resume.

# How to End a Session

---

## Ending a Session

- **After 30 minutes**, the music will slowly fade out.
- **Continue walking for a short period of time after the music stops.** This helps the foot sensors measure your progress.
- Each session will automatically end after 30 minutes (even if a break has been taken).

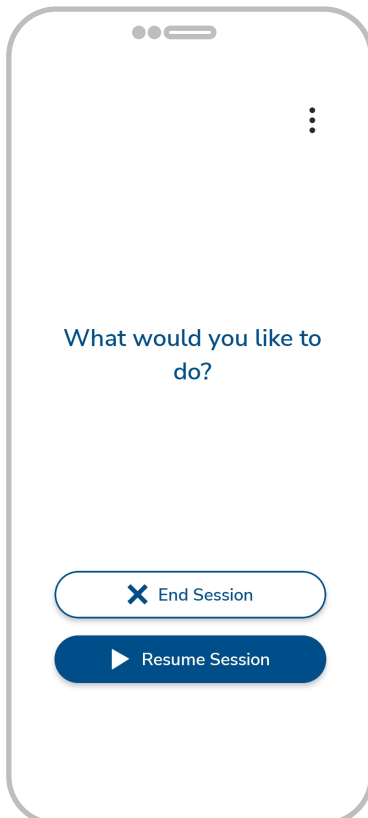If you would like to end a session early, press  at any time and then press


## Remove Foot Sensors

- **Remove the foot sensors from a stable position, such as sitting down.**
- The foot sensors **do not need to be turned off.** They will go to sleep when not in use.
- For your next session, the foot sensors ‘turn on’ just by picking them up. You don’t have to press a button.

## Charge the Move System

- **Keep the system plugged in and charging when you are not using it,** including the foot sensors, touchscreen device, and headset.

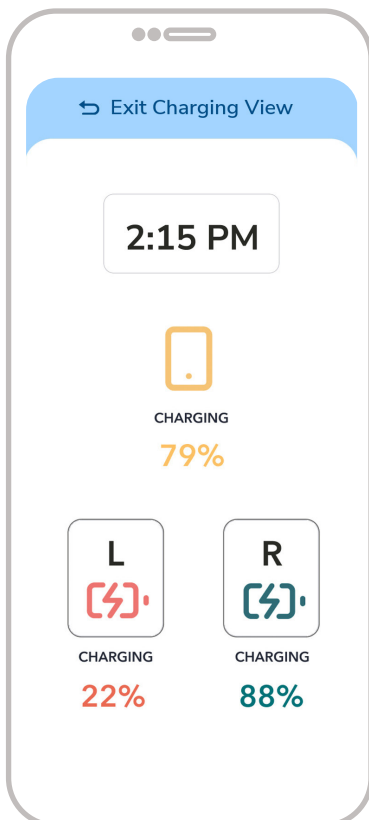During charging, Move automatically uploads all of your session data.

Charging the device also helps to download any updates. That way, Move will be ready for your next session.

## How to Store Moveive

---

### Between Uses, Please:

- **Keep the foot sensors, touchscreen device, headset, and charging equipment somewhere dry** when not in use.
- **Save the shipping box** for returning the system at the end of your prescription

## How to Charge Moveive

---

- ▶ Please see page 11 for instructions on how to charge the system.

# Troubleshooting and Support

If you have any issues using Moveive such as unexpected events or changes in performance, please contact Moveive Support.

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© 2024 MedRhythms, Inc. All rights reserved.
MedRhythms and Move are trademarks of MedRhythms, Inc.
MedRhythms, Inc.
© 2024 MedRhythms, Inc. Todos los derechos reservados.
MedRhythms y Move son marcas comerciales de MedRhythms, Inc.